

User Manual

qPCR-VA a program to analyze qPCR-
data

Murtagh, Justin

Content

1. About the program:.....	2
1.1 When to use this program:.....	2
1.2 What you need to do before using this program:.....	2
2. Requirements before Installation:.....	3
3. Install:	4
3.1 Install with install file:	4
3.2 Manual Install.....	4
4. Usage:	5
5.How to fill out the mask for one sheet/multiple sheets	6
6. Resulting Data	14
7. Example Output	15
8. ERROR HELP	18
8.1 Installation Errors:.....	18
8.2 Double clicking does not work:.....	19
8.3 Working with multiple JREs:	20
9. Example:.....	21

1. About the program:

1.1 When to use this program:

If you ran a qPCR and now want to visualize your results in a plot and calculate the mean efficiency, standard deviation and significance of your data you should use this program. Multiple sheets (resulting from experiments ranging over multiple qPCRs) will be combined. With this program you can compare as many experiment conditions and groups as you like.

The quantification is realized using the $\Delta\Delta\text{Ct}$ -method.

1.2 What you need to do before using this program:

This program is intended to be used as an analytical tool for qPCR-data. The required input is an Excel sheet. You will have to use LinReg on your data before using it with this program, to get the expression levels according to your individual probes. Please remark that you will have to use the group function of LinReg after calling the baselines in order for the program to work properly.

2. Requirements before Installation:

- OS: Windows 64-Bit
- R version 3.4.2 or higher
- (Get newest version at: <https://www.r-project.org/>)
- JRE version 1.8 or higher (Please make sure to only have one JRE installed, or the program will only start from command line!!!)
 - ➔ If you insist to have multiple JRE versions on your computer please refer to Chapter 8.3: Working with multiple JREs (Only recommended if you know what you are doing)

(Get newest version at:

<https://www.oracle.com/technetwork/java/javase/downloads/index.html>)

Optionally

- download R-Studio at:
<https://www.rstudio.com/products/rstudio/download/>
(Only if you intend to change the R-script working in the background)

3. Install:

3.1 Install with install file:

- Right click the file Install.bat and execute it as administrator
- Follow the instructions on the screen
- Make sure all R-packages are installed correctly by starting R and typing
 - `library("rmngb")`
 - `library("reshape2")`
 - `library("ddCt")`
 - `library("ggplot2")`

This step is necessary due to known problems with the automated Installation of R-packages. If Errors occur see the manual installation guide below to install the missing packages.

- qPCR-VA should be ready to start

3.2 Manual Install

- Make sure that you have write permission by starting R as administrator. If this is not possible you need to copy all downloaded packages from your external library (should be in: "C:\Users\Documents\R\win-lib\R3.4\") into the main library of R (should be in: "C:\Program Files\R\R-3.4.2\library")
- Install JRI and rJava by opening R or R-Studio and typing `"install.packages('rJava')"`
- Install rmngb with `"install.packages('rmngb')"`
- Install reshape2 with `"install.packages('reshape2')"`
- Install ddCt with `"source('https://bioconductor.org/biocLite.R')"`
 - `"biocLite('ddCt')"`
- Go to the jri\bin directory -----should be somewhere like: "C:\Program Files\R\R-3.4.2\library\rJava\jri\x64"
 - **!!! Make sure to use the jri.dll from the x64 directory!!!**
- Copy the jri.dll to the jre\bin -----should be somewhere like: "C:\Program Files\Java\jre1.8.0_121\bin"
- Go to the R\bin\x64 directory -----should be somewhere like: "C:\Program Files\R\R-3.4.2\bin\x64"

- Copy all .dll files to the jre\bin -----should be somewhere like "C:\Program Files\Java\jre1.8.0_121\bin"
- Run qPCR-VA.jar
 - Due to known issues with multiple JREs on the computer the JAR-file might not start with a double click. If so, please start the program via the command line using: "java -jar qPCR-VA.jar "
 - ➔ If you are not familiar with the Windows command line see Chapter 8.2: Double clicking does not work

4. Usage:

This program is intended to work with data generated by LinReg. The Resulting Excel sheet should contain 4 Header rows with general Information and then the Data.

Please make sure the Columns are in the following order:

name	indiv_PCR_eff	Amplicon	threshold	mean_PCR_eff	cq	N0

The names must have the format:

Wellplate_Gene_Condition_Group (e.g. A1_TBP_0Gy_15min)

After controlling your Excel sheet start the program and fill in the mask according to the data of your Excel sheet

If you used spacebars or other separators (apart from "_") in the Excel sheet also use them while filling in the Mask

When you filled out the whole mask press 'Apply'. The program will now create a bar plot and show the corresponding data in a Table underneath.

If you would like to change the settings of the plot press 'Change Plot'

- ➔ A dialog will pop up, asking you to change the settings of the plot. Press 'OK' to apply the new Settings.


If you wish to save the Plot and Table press the save button and select the directory where the files should be saved. You will need to select another directory then the one where your excel-sheet is located or the files will be recognized as temporary files and will be deleted!

5.How to fill out the masks

Gene:	<input type="text"/>
Label of the X-axis:	<input type="text"/>
Label of the Y-axis:	<input type="text"/>
Housekeeping Gene:	<input type="text"/>
Control:	<input type="text"/>
Experiment Conditions:	<input type="text"/>
Experiment Groups:	<input type="text"/>


You need to fill in the mask with your data. Attention when filling in the marked sections, for these must be filled exactly how the data is written in Excel. I.e. if you used separators for the gene names or conditions you will have to use the same ones to fill in the mask.

Now choose the type of plot and the resolution you would like to get.

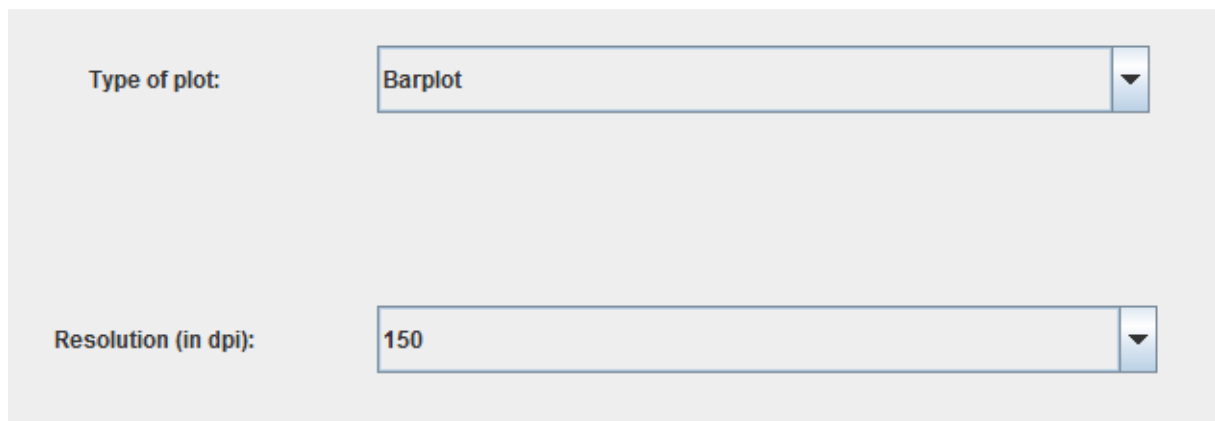

The screenshot shows a configuration window with two dropdown menus. The first dropdown, labeled 'Type of plot:', has 'Barplot' selected. The second dropdown, labeled 'Resolution (in dpi):', has '150' selected.

Now it depends on how many Excel sheets your data is on. (e.g. If you had to run two qPCRs for your Samples the data should be split up on two sheets)
Press the “+”-button until you have enough fields to input your data. (Do not add more sheets than you actually have, or the program will not work)

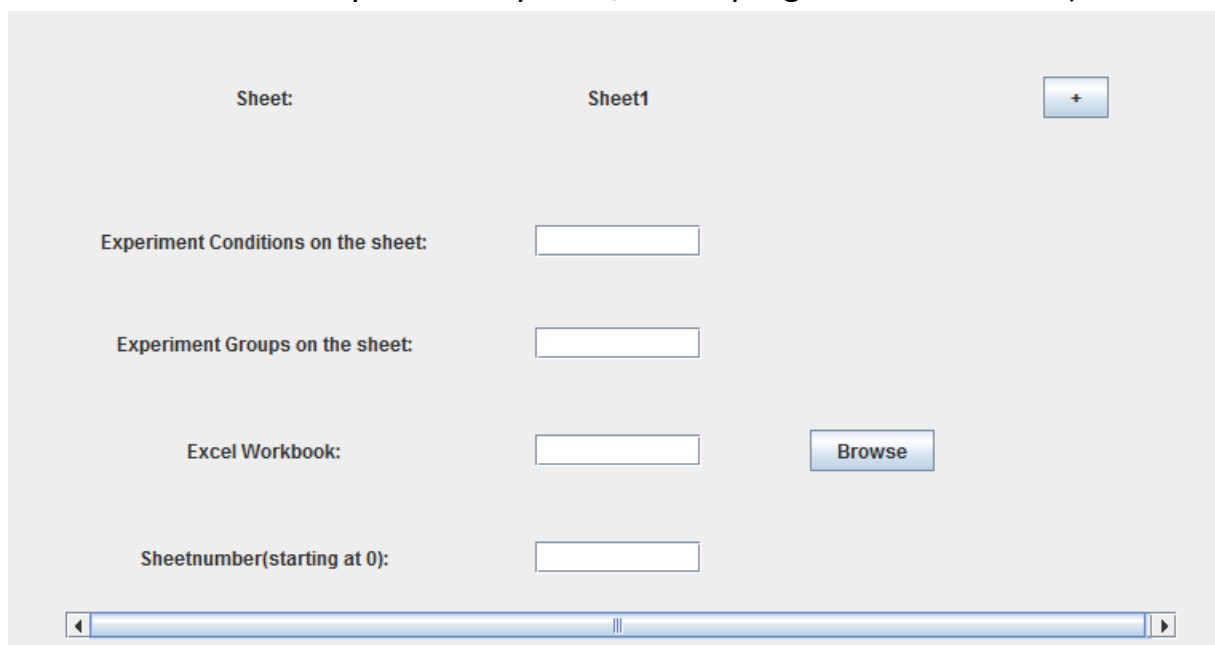

The screenshot shows a configuration window for multiple sheets. At the top, 'Sheet:' is followed by 'Sheet1' and a '+' button. Below this are four input fields: 'Experiment Conditions on the sheet:', 'Experiment Groups on the sheet:', 'Excel Workbook:', and 'Sheetnumber(starting at 0):'. A 'Browse' button is next to the 'Excel Workbook:' field. At the bottom, there is a horizontal scrollbar.

Above is an example for one sheet.


Sheet:	Sheet1		Sheet2	
Experiment Conditions on the sheet:	<input type="text"/>		<input type="text"/>	
Experiment Groups on the sheet:	<input type="text"/>		<input type="text"/>	
Excel Workbook:	<input type="text"/>	<input type="button" value="Browse"/>	<input type="text"/>	<input type="button" value="Browse"/>
Sheetnumber(starting at 0):	<input type="text"/>		<input type="text"/>	

Here an example for two sheets.

Please fill in which data is on which sheet.

If you filled in the whole mask press 'Apply' to start computing.

If you correctly filled out the mask you should get something like this:

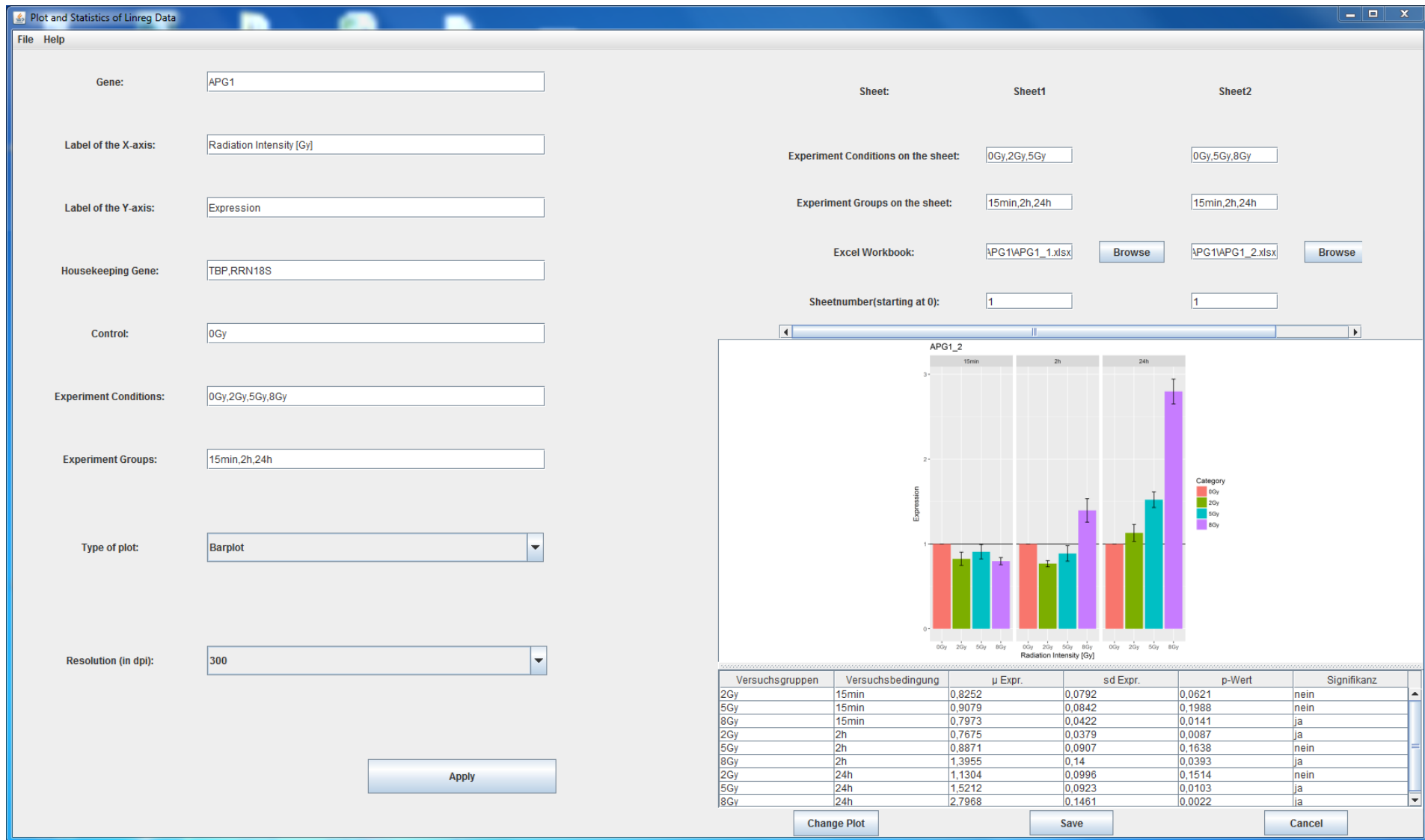

You can now either directly save the generated plot and the table to your desired directory by pressing the save-button, or you can adjust the setting of your plot.

If you wish to adjust your plot press the 'change plot'-button to get the following mask:


Here you can change all possible settings of your plot. Please remark that the scale limit of the Y-Axis needs a minimum and a maximum which must be separated using a “,”. (e.g. 0,4)

If you do not enter a field the default settings will be used.


Click ok to apply the settings to your plot.

!! The plot sometimes does not show right away. If so move the dividing bar between plot and table to refresh the plot.


To perform a statistical Analysis of your data you should do the basic analysis describe above for at least three technical replicates, to guaranty statistical relevant results. After the basic analysis you can then perform the actual statistical analysis with your biological replicates (again at least three samples are recommended for statistical relevance). To do so, go to “File → Statistical Analysis”.

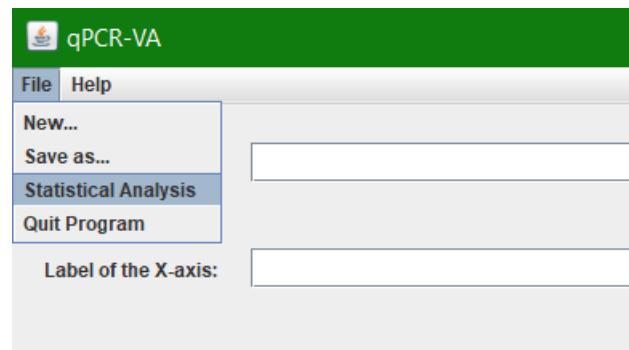

Now choose the files for your analysis using the “browse”- button. (These are the result files from the basic analysis). After selection the files will be shown in the text-field. Please then select whether your data is paired or not.

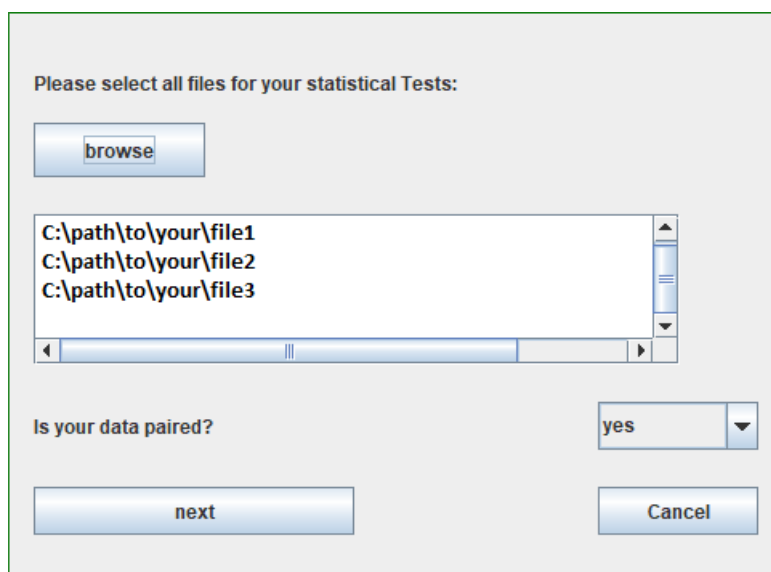

After pressing “next” you are presented a window looking something like this:

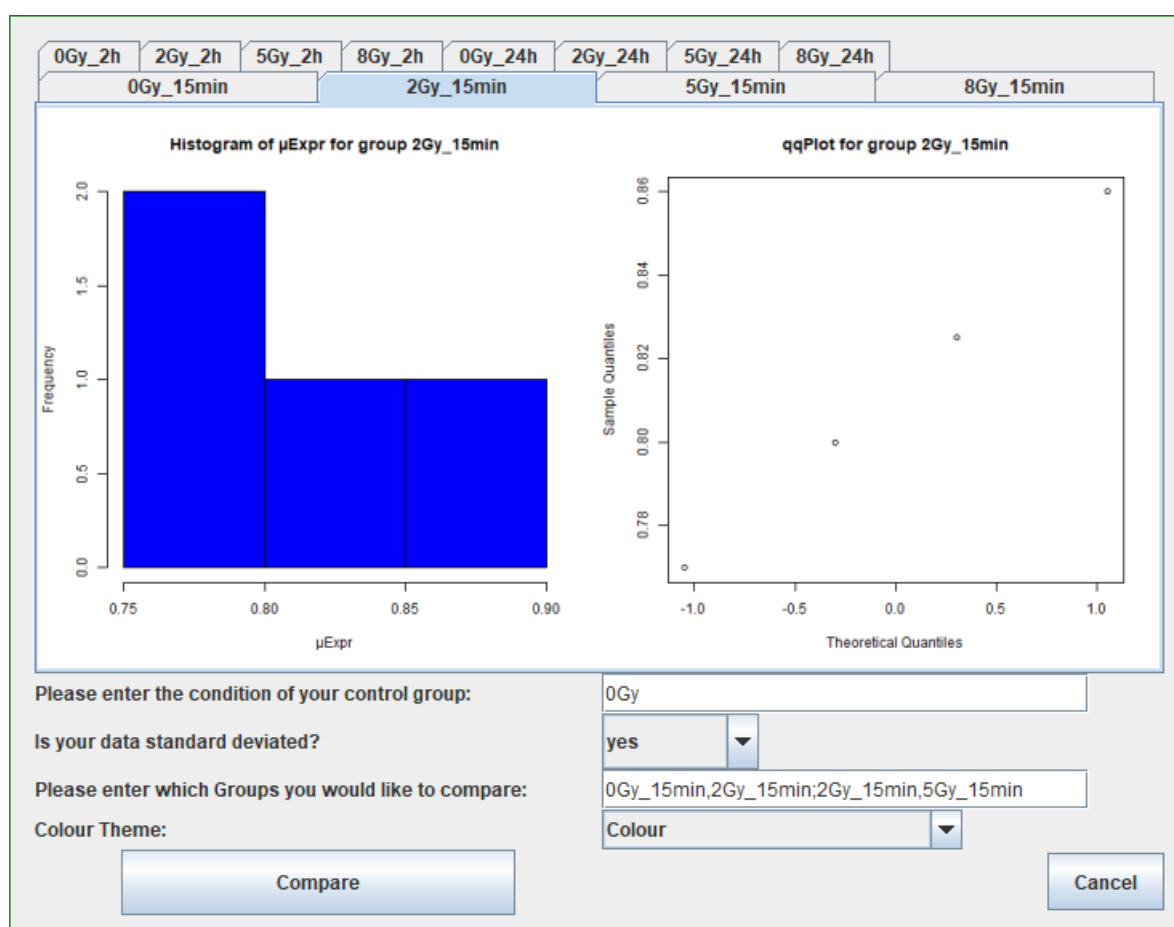


The histograms and qqPlots are presented for each group-condition-pair and are for your assessment on the deviation of your data. If it is standard deviated please select “yes”, if not select no.

To assure the correct statistical test is used in the analysis please enter the control condition which is used for your groups (e.g. 0Gy)

Please also enter which group-condition-pairs you would like to compare. Keep in mind that you can only do one on one comparisons. Separate the pair using a “,” and the different comparisons using “;”. If you want to compare everything with everything you can type “all”. (This is not necessary in most cases and will take quite some time to compute even with small sample numbers).

Finally you can choose whether your resulting plots of the comparison should be in colour or black and white.

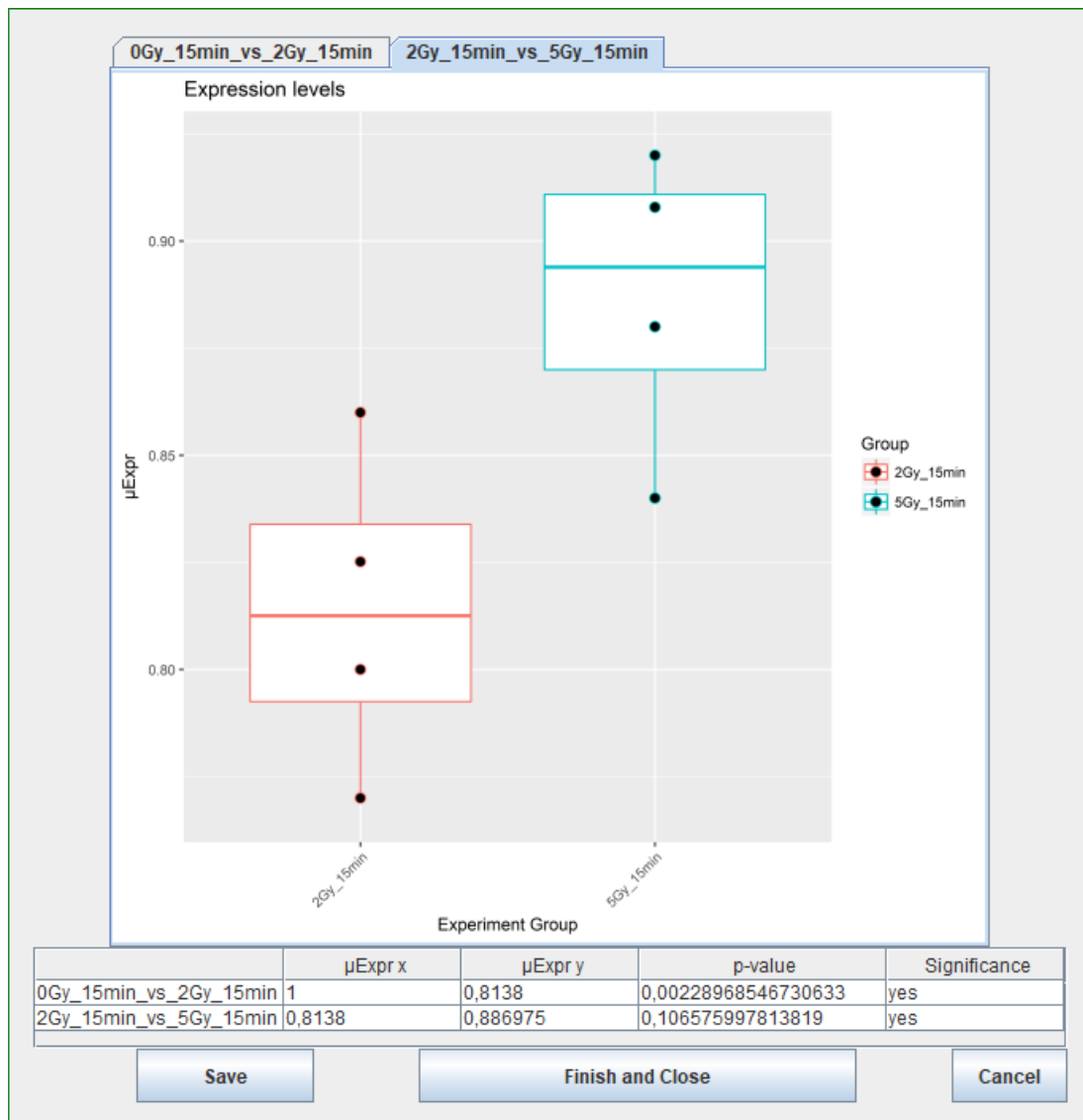

This is the result screen you should now see. The boxplots showing the means and deviations of your data with dots marking each sample. The table below shows the mean expressions of the first and second comparison partner for each comparison as well as the p-value and significance for each comparison.

The statistical test used for computing p-value and significance depends on the answers given:

- Paired and standard deviated: **t-test** (two-sample or one sample against control)
- Paired and not standard deviated: **paired t-test** (two-sample or one sample against control)
- Not paired and standard deviated: **Wilcoxon signed rank test**
- Not paired and not standard deviated: **Wilcoxon rank sum test** (equivalent to the **Mann-Whitney test**)


6. Resulting Data

If everything went correctly you should have two files in your chosen save directory.

1. A '.png' file containing the plot
2. A '.csv' file containing the results of the statistical evaluation of the data


The main reason to use this program should be to find out if the difference in expression levels during the PCR is significant or not, and to gain a visualization of the results.

7. Example Output

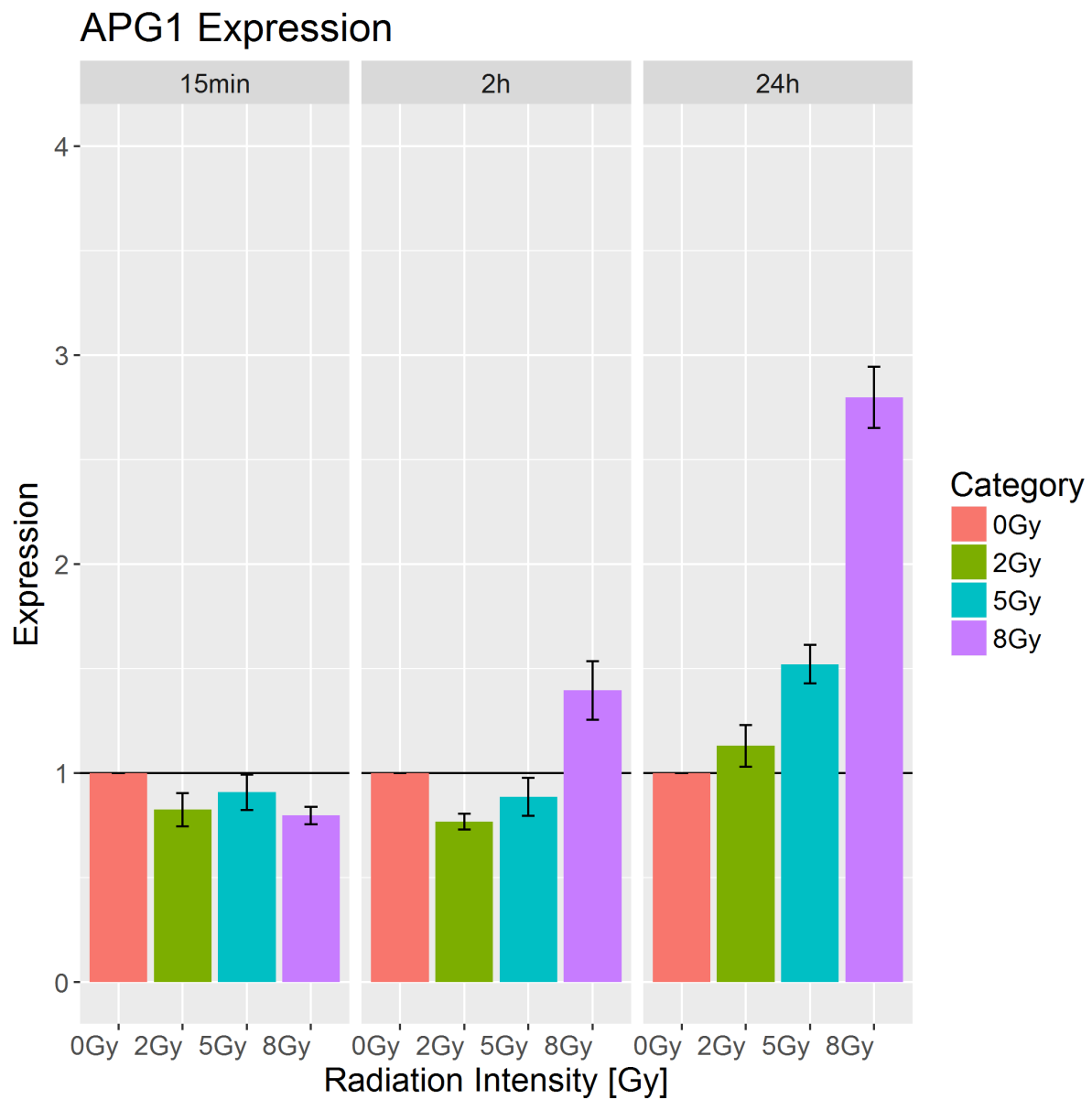

Example with:

- Title = APG1 Expression
- X-label = Radiation Intensity
- Y-label = Expression
- Font size = 15
- Y-limits = 0, 4
- Orientation = Horizontal
- Plot type = Bar plot
- Resolution = 300 dpi

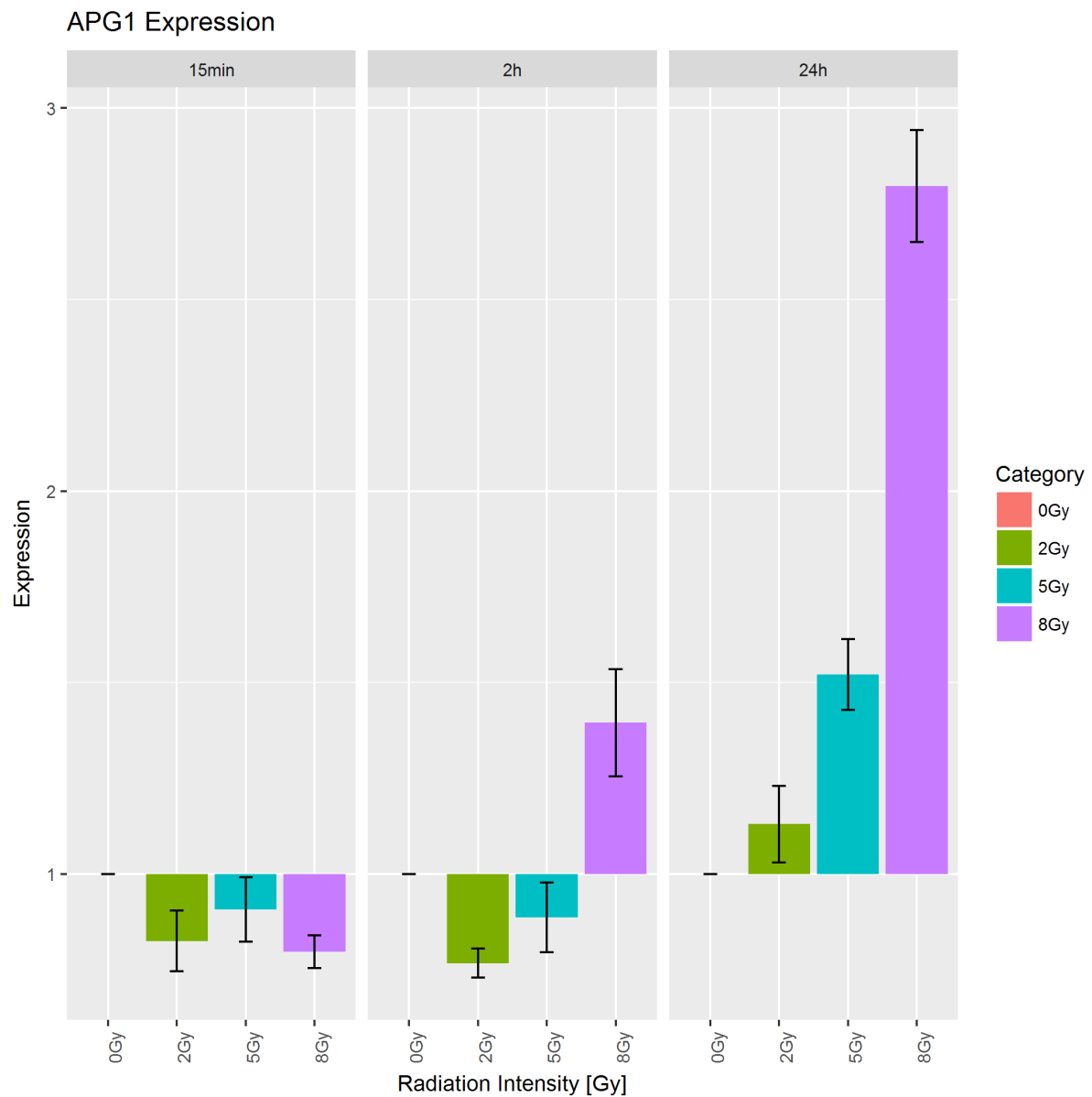

Example with:

- Title = APG1 Expression
- X-label = Radiation Intensity
- Y-label = Expression
- Font size = 15
- Y-limits = 0, 3
- Orientation = Vertical
- Plot type = Bar plot – deviation from mean
- Resolution = 300 dpi

Versuchsgruppen	Versuchsbedingung	μ Expr.	sd Expr.	p-Wert	Signifikanz
2Gy	15min	0,8252	0,0792	0,0621	nein
5Gy	15min	0,9079	0,0842	0,1988	nein
8Gy	15min	0,7973	0,0422	0,0141	ja
2Gy	2h	0,7675	0,0379	0,0087	ja
5Gy	2h	0,8871	0,0907	0,1638	nein
8Gy	2h	1,3955	0,14	0,0393	ja
2Gy	24h	1,1304	0,0996	0,1514	nein
5Gy	24h	1,5212	0,0923	0,0103	ja
8Gy	24h	2,7968	0,1461	0,0022	ja

Example with:

- Two sheets
- Three groups (15min, 2h, 24h)
- Four Conditions (0Gy, 2Gy, 5Gy, 8Gy)
- Statistical test = t-test

8. ERROR HELP

8.1 Installation Errors:

- If you get an Error saying something like:
"Could not move temporary file from 'C:\\\\Path' to 'C:\\\\Path'"

Try typing: `debug(utils::unpackPkgZip)` and then installing the package again

➔ You will have to press Enter a lot until it should say:
"Package downloaded successfully ...something... check sum MD5"

Sometimes it can be necessary to repeat the step mentioned above to get rid of any issues


- Check if all necessary packages are downloaded correctly, by typing:

```
library("rmngb")  
library("reshape2")  
library("ddCt")  
library("ggplot2")
```


➔ If there is a Warning like:
Could not find package "NameOfPackage",
try installing it manually with `"install.packages('NameOfPackage')"`

8.2 Double clicking does not work:

- Open the command line by pressing the windows button + R at the same time. You should see this:

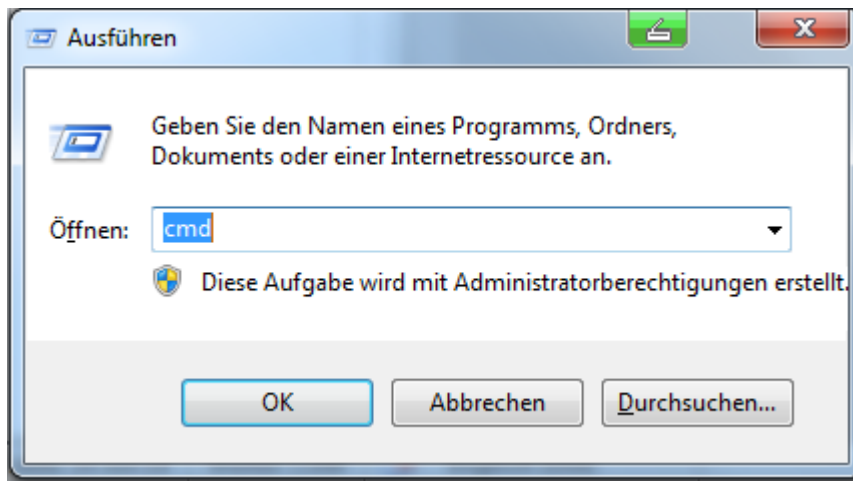

- Type in “cmd” and press Enter
- The command line should come up and show you something like this:

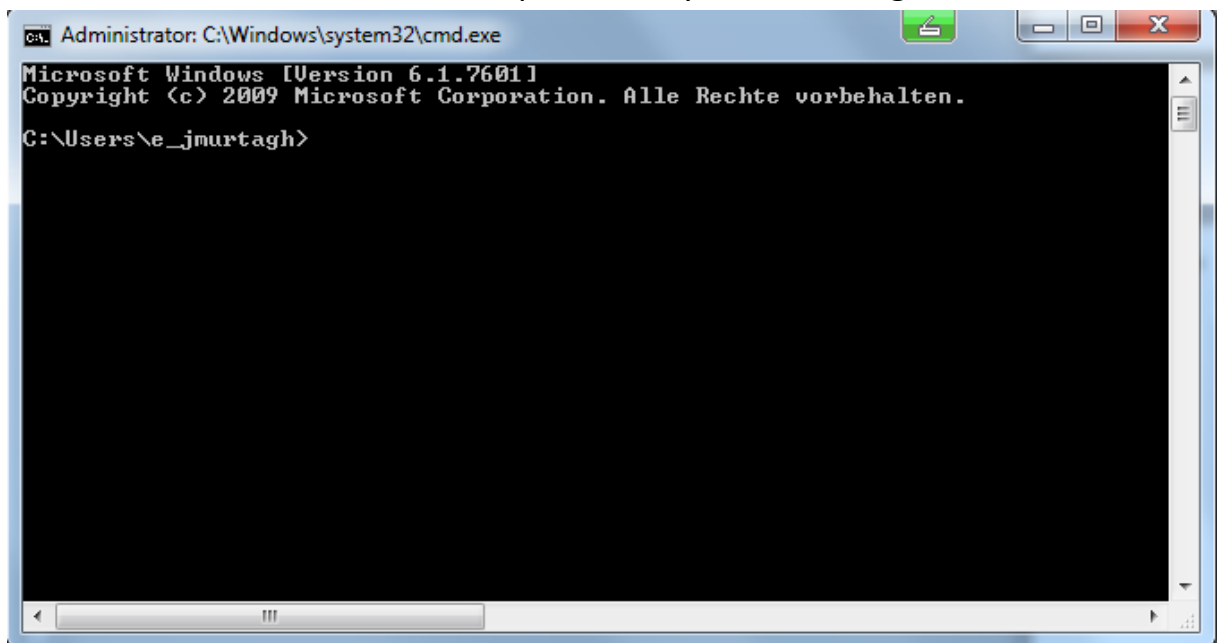

- Now you need to change to the directory of the program
 - E.g. if you installed it in a directory on the desktop type in:
 - `cd C:\Users\YourUser\Desktop\qPCR-VA`
 - (instead of “YourUser” type in the name of the User your working on and instead of “qPCR-VA” the name of the directory the program is in)
- You should now be in the right directory and be able to start the program by typing: “java -jar qPCR-VA.jar”

8.3 Working with multiple JREs:

Change the Registry:

First open the Windows command line by [Windows-button] + [R] and type regedit. Confirm by pressing [Enter].

You now may have to enter the administrator password or you have to confirm the security check of the User Access Control (UAC).

Next the Registry-Editor should open. On the left hand side you will find the code key.

You now need to click until you reach:

HKEY_CLASSES_ROOT\jarfile\shell\open\command

or (depending on the Windows-Version)

HKEY_CLASSES_ROOT\jar_auto_file\shell\open\command

Now open the key editor by double clicking the name standard on the right side of the window.

The value of the key consists of the path to the interpreter and the Wildcards, who act as placeholders for the given data-name:

"C:\Program Files (x86)\Java\jre8\bin\javaw.exe" -jar "%1" %*

The path to the JRE has to be read or copied from the system and then, respective to the template in the line with grey background, pasted to the key-editor. Pressing OK confirms the edit. The Registry-Editor can now be closed.

The Jar-Files should now start by double clicking.

9. Example:

An example Excel-sheet is provided within the program directory.

To test if your installation and setup was successful it is recommended to try out the program with the Example data.

Here is the main input you need to provide:

Gene: APG1

Label of X-axis: Dose

Label of Y-axis: Expression

Housekeeping Gene: TBP,RRN18S

Control: 0Gy

Experiment Conditions: 0Gy,2Gy,5Gy,8Gy

Experiment Groups: 15min,2h,24h

Type of plot: Barplot

Resolution (in dpi): 600

Then open to sheet inputs by pressing the “+”-button two times.

For sheet 1 enter:

Experiment Conditions: 0Gy,2Gy,5Gy

Experiment Groups: 15min,2h,24h

Excel Workbook: Select Example File 1

Sheet number: 0

For sheet2 enter:

Experiment Conditions: 0Gy,5Gy,8Gy

Experiment Groups: 15min,2h,24h

Excel Workbook: Select Example File 2

Sheet number: 0