

Universitätsmedizin Mainz

User Manual

qPCR-VA a program to analyze qPCR-
data

Murtagh, Justin

2017

Content

1. About the program:.....	2
1.1 When to use this program:.....	2
1.2 What you need to do before using this program:.....	2
2. Requirements before Installation:.....	3
3. Install:	4
3.1 Install with install file:	4
3.2 Manual Install.....	4
4. Usage:	6
5.How to fill out the mask for one sheet/multiple sheets	7
6. Resulting Data	11
7. Example Output	12
8. ERROR HELP	15
8.1 Installation Errors:.....	15
8.2 Double clicking does not work:.....	16
8.3 Working with multiple JREs:	17
9. Example:.....	18

1. About the program:

1.1 When to use this program:

If you ran a qPCR and now want to visualize your results in a plot and calculate the mean efficiency, standard deviation and significance of your data you should use this program. Multiple sheets (resulting from experiments ranging over multiple qPCRs) will be combined. With this program you can compare as many experiment conditions and groups as you like.

The quantification is realized using the $\Delta\Delta C_t$ -method.

1.2 What you need to do before using this program:

This program is intended to be used as an analytical tool for qPCR-data. The required input is an Excel sheet. You will have to use LinReg on your data before using it with this program, to get the expression levels according to your individual probes. Please remark that you will have to use the group function of LinReg after calling the baselines in order for the program to work properly.

2. Requirements before Installation:

- OS: Windows 64-Bit
- R version 3.4.2 or higher
- (Get newest version at: <https://www.r-project.org/>)
- JRE version 1.8 or higher (Please make sure to only have one JRE installed, or the program will only start from command line!!!)
 - ➔ If you insist to have multiple JRE versions on your computer please refer to Chapter 8.3: Working with multiple JREs (Only recommended if you know what you are doing)

(Get newest version at:

<https://www.oracle.com/technetwork/java/javase/downloads/index.html>)

Optionally

- download R-Studio at:
<https://www.rstudio.com/products/rstudio/download/>
(Only if you intend to change the R-script working in the background)

3. Install:

3.1 Install with install file:

- Start R and type: `source("C:\\Path\\to\\Inst.R")`
 - Make sure that you have write permission by starting R as administrator. If this is not possible you need to copy all downloaded packages from your external library (should be in: `"C:\\Users\\Documents\\R\\win-lib\\R3.4\\"`) into the main library of R (should be in: `"C:\\Program Files\\R\\R-3.4.2\\library"`)
- Double click the file Install.bat
- Follow the instructions on the screen
- PlotAndStats should be ready to start

3.2 Manual Install

- Install JRI and rJava by opening R or R-Studio and typing `"install.packages('rJava')"`
- Install rmngb with `"install.packages('rmngb')"`
- Install reshape2 with `"install.packages('reshape2')"`
- Install ddCt with `"source('https://bioconductor.org/biocLite.R')"`
 - `"biocLite('ddCt')"`
- Go to the `jri\\bin` directory -----should be somewhere like:
`"C:\\Program Files\\R\\R-3.4.2\\library\\rJava\\jri\\x64"`
 - **!!! Make sure to use the jri.dll from the x64 directory!!!**
- Copy the `jri.dll` to the `jre\\bin` -----should be somewhere like:
`"C:\\Program Files\\Java\\jre1.8.0_121\\bin"`
- Go to the `R\\bin\\x64` directory -----should be somewhere like:
`"C:\\Program Files\\R\\R-3.4.2\\bin\\x64"`
- Copy all `.dll` files to the `jre\\bin` -----should be somewhere like
`"C:\\Program Files\\Java\\jre1.8.0_121\\bin"`

- Run PlotAndStats.jar
 - Due to known issues with multiple JREs on the computer the JAR-file might not start with a double click. If so, please start the program via the command line using: “java -jar PlotAndStats.jar “
 - ➔ If you are not familiar with the Windows command line see Chapter 8.2: Double clicking does not work

4. Usage:

This program is intended to work with data generated by LinReg. The Resulting Excel sheet should contain 4 Header rows with general Information and then the Data.

Please make sure the Columns are in the following order:

name	indiv_PCR_eff	Amplicon	threshold	mean_PCR_eff	cq	N0
-------------	----------------------	-----------------	------------------	---------------------	-----------	-----------

The names must have the format:

Wellplate_Gene_Condition_Group (e.g. A1_TBP_0Gy_15min)

After controlling your Excel sheet start the program and fill in the mask according to the data of your Excel sheet

If you used spacebars or other separators (apart from "_") in the Excel sheet also use them while filling in the Mask

When you filled out the whole mask press 'Apply'. The program will now create a bar plot and show the corresponding data in a Table underneath.

If you would like to change the settings of the plot press 'Change Plot'

➔ A dialog will pop up, asking you to change the settings of the plot.
Press 'OK' to apply the new Settings.

If you wish to save the Plot and Table press the save button and select the directory where the files should be saved. You will need to select another directory then the one where your excel-sheet is located or the files will be recognized as temporary files and will be deleted!

5.How to fill out the mask for one sheet/multiple sheets

Gene:	<input type="text"/>
Label of the X-axis:	<input type="text"/>
Label of the Y-axis:	<input type="text"/>
Housekeeping Gene:	<input type="text"/>
Control:	<input type="text"/>
Experiment Conditions:	<input type="text"/>
Experiment Groups:	<input type="text"/>

You need to fill in the mask with your data. Attention when filling in the marked sections, for these have to be filled exactly how the data is written in Excel. I.e. if you used separators for the gene names or conditions you will have to use the same ones to fill in the mask.

Now choose the type of plot and the resolution you would like to get.

Type of plot:	<input type="text" value="Barplot"/>
Resolution (in dpi):	<input type="text" value="150"/>

Now it depends on how many Excel sheets your data is on. (e.g. If you had to run two qPCRs for your Samples the data should be split up on two sheets)
Press the “+”-button until you have enough fields to input your data. (Do not add more sheets than you actually have or the program will not work)

The screenshot shows a software interface for configuring data input for a single sheet. At the top, there is a label 'Sheet:' followed by 'Sheet1' and a blue button with a '+' sign. Below this, there are four rows of input fields: 'Experiment Conditions on the sheet:', 'Experiment Groups on the sheet:', 'Excel Workbook:', and 'Sheetnumber(starting at 0):'. Each row has a corresponding text input box. To the right of the 'Excel Workbook:' input box is a blue 'Browse' button. At the bottom of the interface is a horizontal scrollbar.

Above is an example for one sheet.


The screenshot shows the software interface configured for two sheets. The top section has 'Sheet:' followed by 'Sheet1' and 'Sheet2'. Below this, there are four rows of input fields for each sheet: 'Experiment Conditions on the sheet:', 'Experiment Groups on the sheet:', 'Excel Workbook:', and 'Sheetnumber(starting at 0):'. Each row has two text input boxes, one for each sheet. To the right of the 'Excel Workbook:' input boxes are two blue 'Browse' buttons. At the bottom of the interface is a horizontal scrollbar.

Here an example for two sheets.

Please fill in which data is on which sheet.

If you filled in the whole mask press 'Apply' to start computing.

If you correctly filled out the mask you should get something like this:

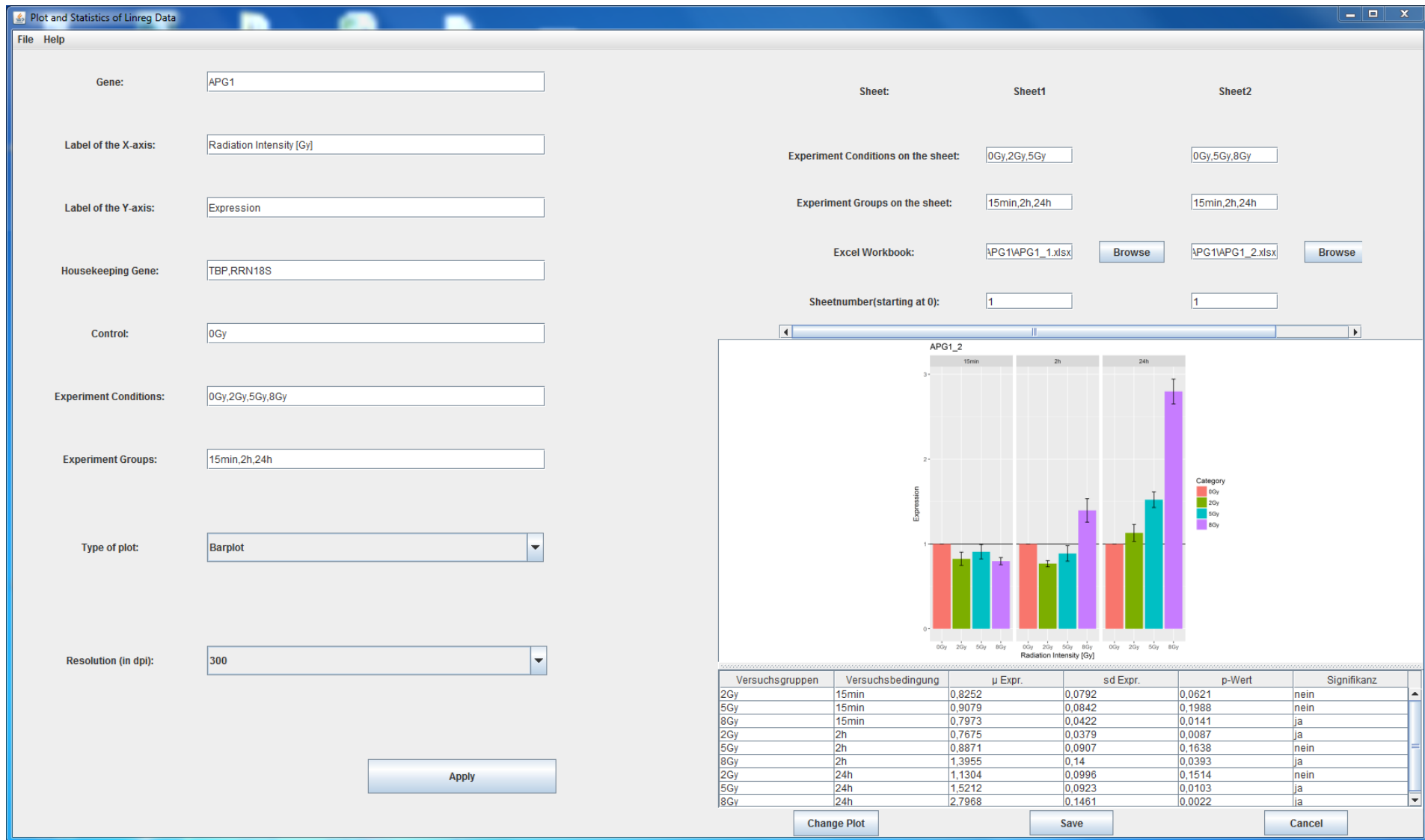

You can now either directly save the generated plot and the table to your desired directory by pressing the save-button, or you can adjust the setting of your plot.

If you wish to adjust your plot press the 'change plot'-button to get the following mask:

Here you can change all possible settings of your plot. Please remark that the scale limits of the Y-Axis needs a minimum and a maximum which has to be separated using a “,”. (e.g. 0,4)

If you do not enter a field the default settings will be used.

Click ok to apply the settings to your plot.

!! The plot sometimes does not show right away. If so move the dividing bar between plot and table to refresh the plot.


6. Resulting Data

If everything went correctly you should have two files in your chosen save directory.

1. A '.png' file containing the plot
2. A '.csv' file containing the results of the statistical evaluation of the data


The main reason to use this program should be to find out if the difference in expression levels during the PCR is significant or not, and to gain a visualization of the results.

7. Example Output

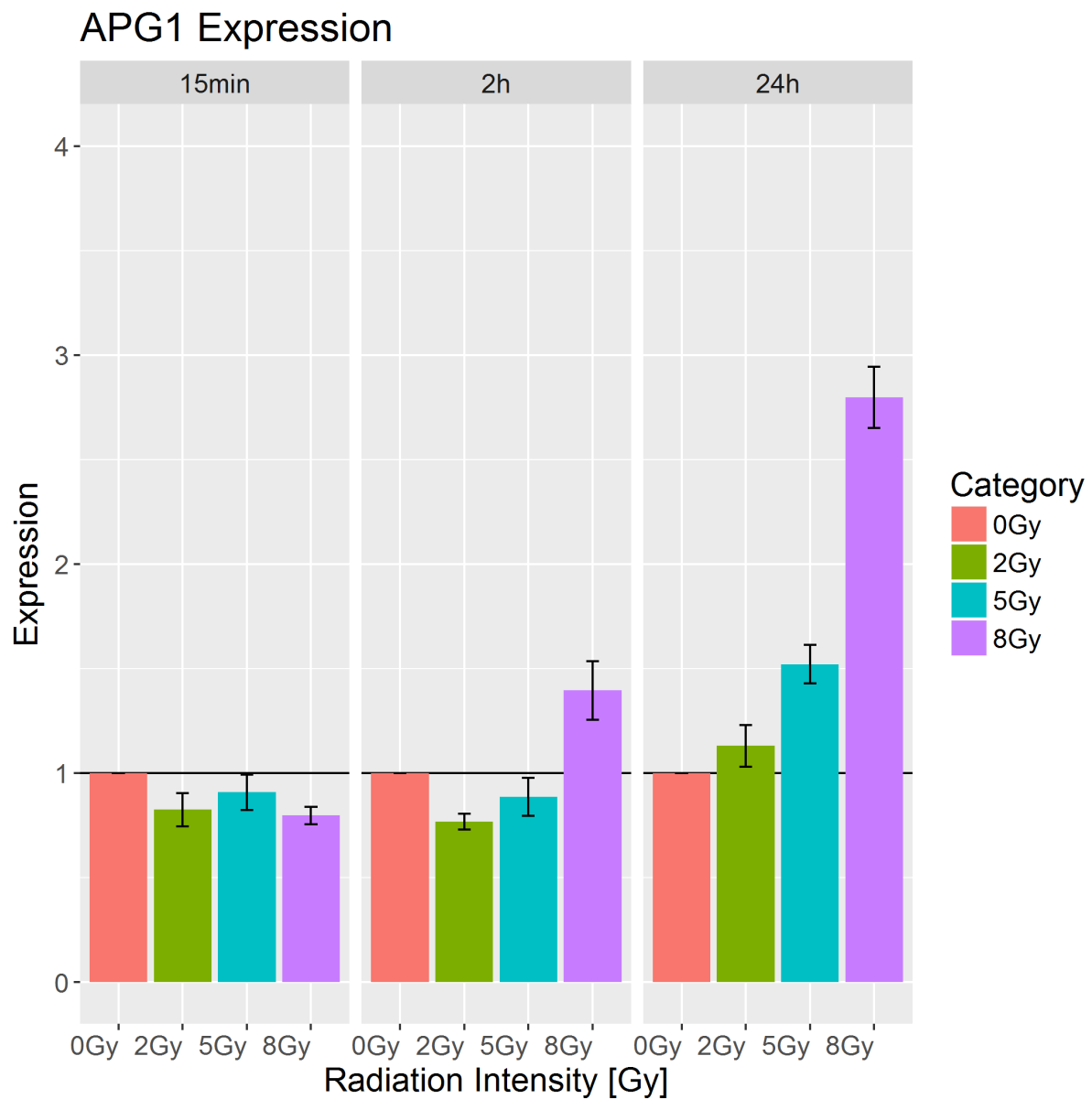

Example with:

- Title = APG1 Expression
- X-label = Radiation Intensity
- Y-label = Expression
- Font size = 15
- Y-limits = 0, 4
- Orientation = Horizontal
- Plot type = Bar plot
- Resolution = 300 dpi

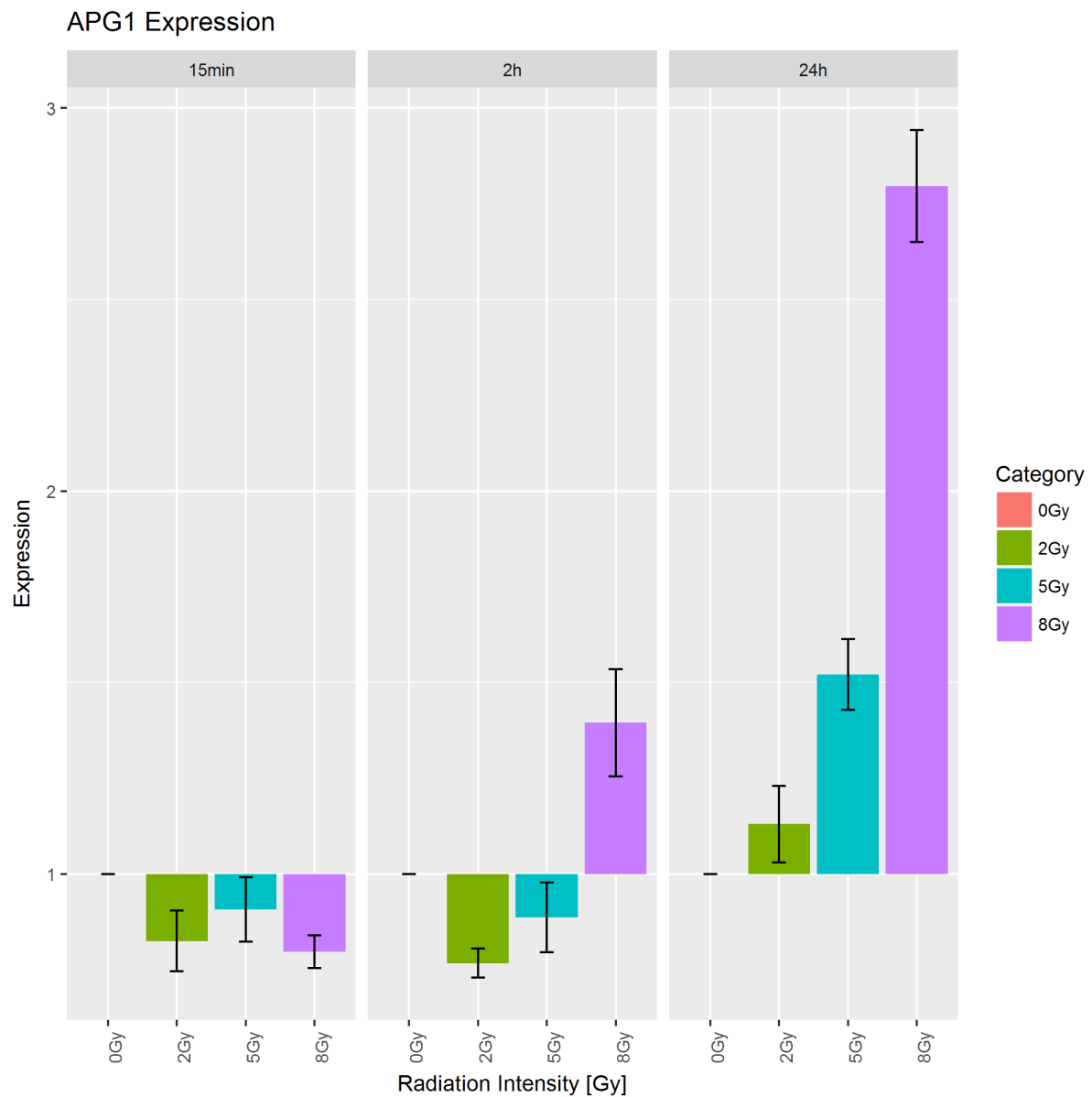

Example with:

- Title = APG1 Expression
- X-label = Radiation Intensity
- Y-label = Expression
- Font size = 15
- Y-limits = 0, 3
- Orientation = Vertical
- Plot type = Bar plot – deviation from mean
- Resolution = 300 dpi

Versuchsgruppen	Versuchsbedingung	μ Expr.	sd Expr.	p-Wert	Signifikanz
2Gy	15min	0,8252	0,0792	0,0621	nein
5Gy	15min	0,9079	0,0842	0,1988	nein
8Gy	15min	0,7973	0,0422	0,0141	ja
2Gy	2h	0,7675	0,0379	0,0087	ja
5Gy	2h	0,8871	0,0907	0,1638	nein
8Gy	2h	1,3955	0,14	0,0393	ja
2Gy	24h	1,1304	0,0996	0,1514	nein
5Gy	24h	1,5212	0,0923	0,0103	ja
8Gy	24h	2,7968	0,1461	0,0022	ja

Example with:

- Two sheets
- Three groups (15min, 2h, 24h)
- Four Conditions (0Gy, 2Gy, 5Gy, 8Gy)
- Statistical test = t-test

8. ERROR HELP

8.1 Installation Errors:

- If you get an Error saying something like:
"Could not move temporary file from 'C:\\\\Path' to 'C:\\\\Path'"

Try typing: `debug(utils::unpackPkgZip)` and then installing the package again

➔ You will have to press Enter a lot until it should say:
"Package downloaded successfully ...something... check sum MD5"

Sometimes it can be necessary to repeat the step mentioned above to get rid of any issues


- Check if all necessary packages are downloaded correctly, by typing:

```
library("rmnbg")  
library("reshape2")  
library("ddCt")  
library("ggplot2")
```


➔ If there is a Warning like:
Could not find package "NameOfPackage",
try installing it manually with `"install.packages('NameOfPackage')"`

8.2 Double clicking does not work:

- Open the command line by pressing the windows button + R at the same time. You should see this:

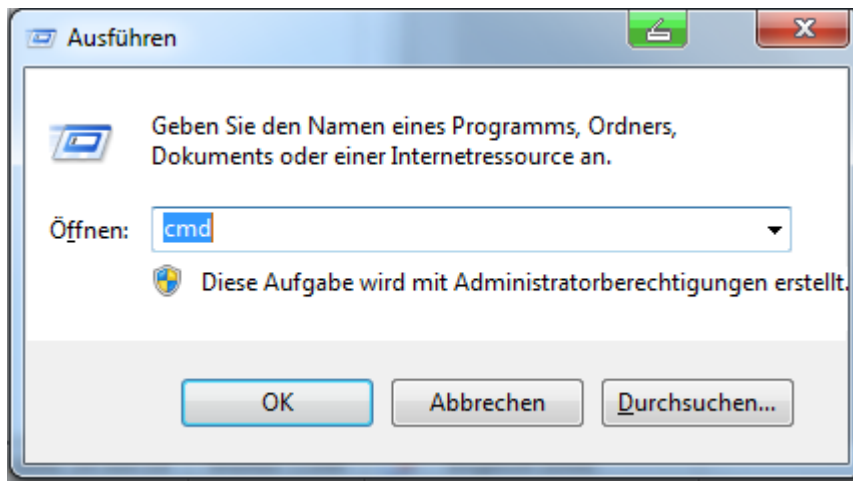

- Type in “cmd” and press Enter
- The command line should come up and show you something like this:

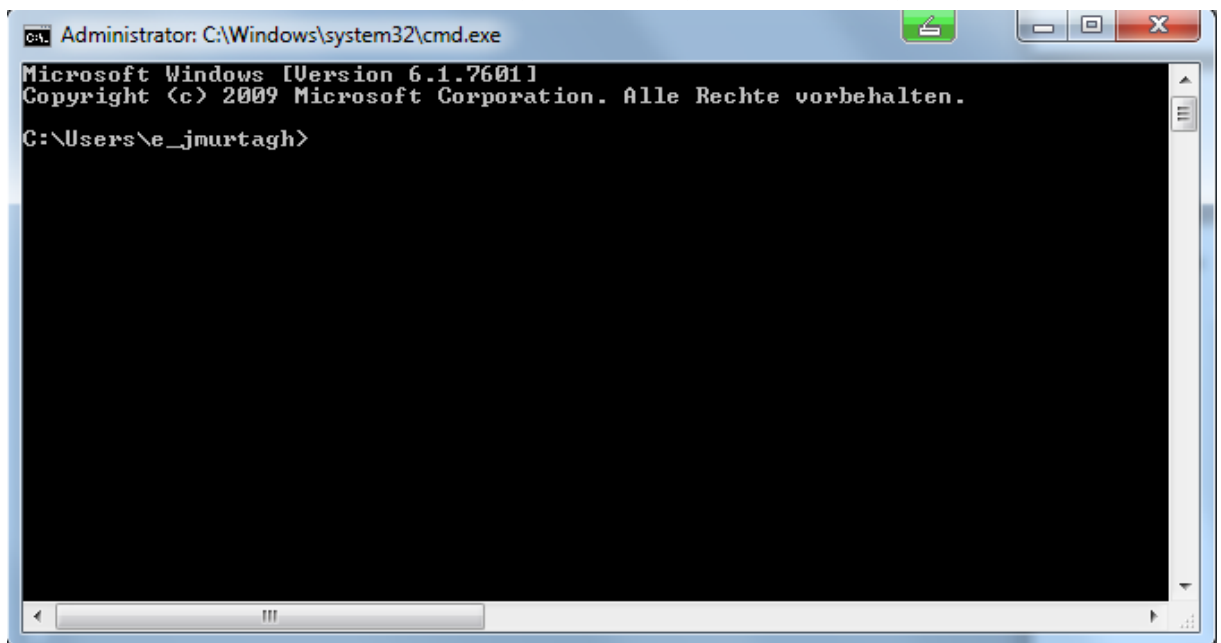

- Now you need to change to the directory of the program
 - E.g. if you installed it in a directory on the desktop type in:
 - `cd C:\Users\YourUser\Desktop\PlotAndStats`
 - (instead of “YourUser” type in the name of the User your working on and instead of “PlotAndStats” the name of the directory the program is in)
- You should now be in the right directory and be able to start the program by typing: “java -jar PlotAndStats.jar”

8.3 Working with multiple JREs:

Change the Registry:

First open the Windows command line by [Windows-button] + [R] and type regedit. Confirm by pressing [Enter].

You now may have to enter the administrator password or you have to confirm the security check of the User Access Control (UAC).

Next the Registry-Editor should open. On the left hand side you will find the code key.

You now need to click until you reach:

HKEY_CLASSES_ROOT\jarfile\shell\open\command

or (depending on the Windows-Version)

HKEY_CLASSES_ROOT\jar_auto_file\shell\open\command

Now open the key editor by double clicking the name standard on the right side of the window.

The value of the key consists of the path to the interpreter and the Wildcards, who act as placeholders for the given data-name:

"C:\Program Files (x86)\Java\jre8\bin\javaw.exe" -jar "%1" %*

The path to the JRE has to be read or copied from the system and then, respective to the template in the line with grey background, pasted to the key-editor. Pressing OK confirms the edit. The Registry-Editor can now be closed.

The Jar-Files should now start by double clicking.

9. Example:

An example Excel-sheet is provided within the program directory.

To test if your installation and setup was successful it is recommended to try out the program with the Example data.

Here is the main input you need to provide:

Gene: APG1

Label of X-axis: Dose

Label of Y-axis: Expression

Housekeeping Gene: TBP,RRN18S

Control: 0Gy

Experiment Conditions: 0Gy,2Gy,5Gy,8Gy

Experiment Groups: 15min,2h,24h

Type of plot: Barplot

Resolution (in dpi): 600

Then open to sheet inputs by pressing the “+”-button two times.

For sheet 1 enter:

Experiment Conditions: 0Gy,2Gy,5Gy

Experiment Groups: 15min,2h,24h

Excel Workbook: Select the Example File

Sheet number: 0

For sheet2 enter:

Experiment Conditions: 0Gy,5Gy,8Gy

Experiment Groups: 15min,2h,24h

Excel Workbook: Select the Example File

Sheet number: 1